

Bluetooth Solution Kit

Wireless sound transmission directly to your hearing instrument from both stationary and mobile phone

Wireless sound transmission

Bluetooth Solution Kit facilitates mobile phone communication by wireless sound transmission to the hearing instrument by aid of the accompanying neck loop.

If you have a built-in MP3 player in your mobile phone you can use that function in a corresponding way. You can then switch between talking in the mobile phone and e.g. listening to music. By using the T mode in your hearing instrument you minimise background noise, and it becomes considerably easier to hear!

What equipment is needed?

- A hearing instrument with built-in telecoil, which the very most hearing instruments have.
- Bluetooth in your mobile phone, which most modern mobile phones have.

Functions in the Bluetooth device

- Buttons with clear symbols on the Bluetooth device
 - Answering calls/ending calls
 - Volume up/down
 - For music: play/pause/next track/previous track
- Accompanying charger
- Up to 8 hours talking time
- Up to 7 hours music playing
- Up to 230 hours standby time
- Splash proof
- Possibility to connect two units simultaneously

Accessories for fixed/IP telephony

If you would like to use Bluetooth Solution Kit also for fixed phones you can buy the adapter "Jabra A7010 Bluetooth Hub". It is suited for all fixed phone networks and different phone models.

If you have IP telephony on your computer, you can add on the accessory "Jabra A330 BT Dongle", a Bluetooth dongle that is inserted in the USB port on your computer. Note the great advantage that you can connect Bluetooth Solution Kit to both a mobile phone and a fixed/IP phone in parallel!



Neck loop, 2-channel

This product is CE marked according to the EU Medical Device Directive 93/42/EEC.

Bluetooth unit (BT3030)

This product is CE marked according to the provisions of the R & TTE directive (99/5/EC).



Ordering numbers

Bluetooth Solution Kit	2002205
Jabra A7010 Bluetooth hub	2002208
Jabra A330 Bluetooth Dongle	2002207
Hook Switch GN 1000 RHL	1110004



Bluetooth Solution Kit



Hub A7010



Dongle A330



Hook Switch GN1000 RHL

9950113 Issue 2, 2010-01 We reserve the right to make changes in performance and specification without prior notice.

MANUFACTURER

in the sense of EU Medical Device Directive 93/42/EEC:



GN ReSound A/S

Lautrupbjerg 7
DK-2750 Ballerup
Denmark



Any issues relating to the EU Medical Device Directive 93/42/EEC should be directed to GN ReSound A/S.

Service and product-related questions: contact GN ReSound AB ALD Division.

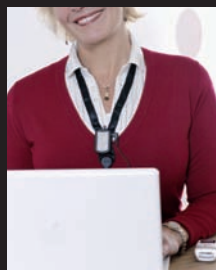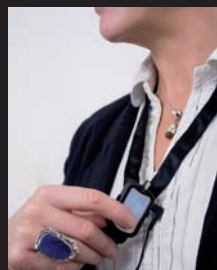

www.gnresound-ald.com

GN ReSound AB ALD Division
Box 12874
SE-112 98 STOCKHOLM
Sweden

Tel. 08-441 16 70
Fax 08-441 16 80
info.ald@gnresound.se

ReSound

rediscover hearing